
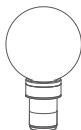
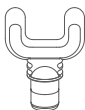


KiCA


Vibration Percussion Device

Instructions V2.1 Model: K2

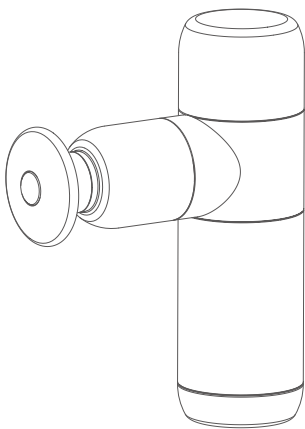

Overview

Main Body

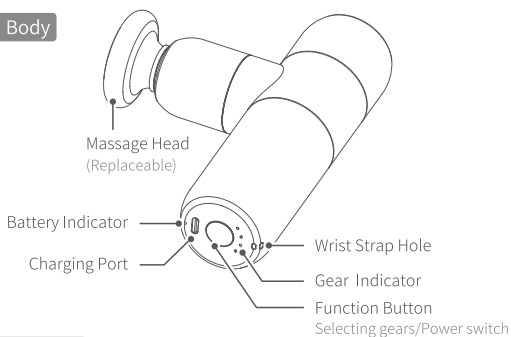

Massage Head

Magnetic



Flat Head

Massage flat major muscle groups on the back

Magnetic



Cone Head

Percuss deep tissues like meridians and acupoint, etc.



Spherical Head

Relax soft tissues and sensitive muscles




U-shaped Head

Massage neck and tendon


Start to Use

Install/ Replace Massage Head

Install a massage head into the hole in front of the device according to your need.

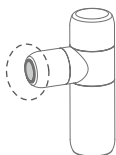


Installation Rendering:

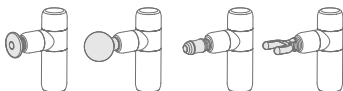

Power ON

Long press the function button.

! Do not use it while charging.

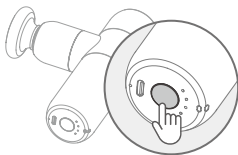

Power OFF

Single tap function button to switch to standby mode or long press function button.

Selecting Gears

Single tap function button to select gears after powering on.

The product will automatically power off after enter standby mode for 2 seconds.

Gears	Rated Rotate Speed (rpm)
First gear	1200
Second gear	2000
Third gear	3200
Intelligent massage	1200~3200
Standby mode	0

Intelligent Massage

This product has intelligent massage function. Under this gear the product will automatically give a circular massage between high and low speed. The highest speed is 3200 rpm and the lowest speed is 1200 rpm.

Intelligent Protection

This product has intelligent protection function. In order to prevent muscle injury caused by massaging for a long time, gear indicator will turn off (battery indicator remain) after working for 10 minutes continuously, and it can be restarted by pressing the function button. If there is no operation within 5 seconds, it will power off automatically.

Locked Motor Protection

This product has locked motor protection function. If the rotating is blocked for 3 seconds, the product will enter standby mode automatically. Users can restart it by pressing the function button. If no operation in 10 seconds, the product will power off automatically.

Indicator


Battery Indicator

Light	Status	Description
Red	████████	Charging
Green	████████	Charge complete
Green	████████	Fully charged
Green	███	Medium level of battery
Green	●●●●●	Low power, need recharging

Gear Indicator

Light	Status	Description
Blue	●○○○	First gear
Blue	●●○○	Second gear
Blue	●●●○	Third gear
Blue	●●●●	Intelligent massage

Charging

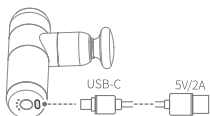

Do not use it while charging

Specifications

Product name:

KiCA Vibration Percussion Device

Model: K2

Operating Voltage Range: 10.6V-12.6V

Battery Capacity: 2000mAh

Charging Time: >2.5 Hours


Weight (Reference value): 485g (Main Body)

Size: 90.5 x 144.6 x 44mm

Theoretical Battery Life: About 11 Hours

List

- | | | | |
|-----------------|----|---------------|----|
| ① Main Body | x1 | ⑥ USB-C Cable | x1 |
| ② Flat Head | x1 | ⑦ Storage Bag | x1 |
| ③ Cone Head | x1 | ⑧ Wrist Strap | x1 |
| ④ Spherica Head | x1 | ⑨ Manual | x1 |
| ⑤ U-shaped Head | x1 | | |

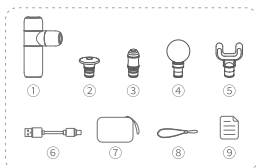

! Contraindications

- (1) Cancers.
- (2) Metal bone anastomosis (do not use the device directly in the fracture area).
- (3) Osteoporosis.
- (4) Vein thrombosis and pulsation.
- (5) Hypertension.
- (6) Pregnancy (do not use the device in close proximity to the abdomen).
- (7) Kidney stones (do not use the device near the kidneys).
- (8) Varicose veins.
- (9) Acute inflammation.
- (10) Do not massage one episode for more than 4-5 minutes.
- (11) Do not massage the bony area.



Do not press the massage head firmly to block the motor rotation while using, otherwise it will reduce the service life of the product!

Disclaimer

Thanks for using the product. The information in this document affects your safety, legal rights, and responsibilities. Read the entire document carefully to ensure proper configuration before using. Failure to read and follow instructions and warnings in this document may result in serious injury to yourself or others, or damage your products or other objects in the vicinity.

By using this product, you hereby signify that you have read this disclaimer and warning carefully and that you understand and agree to abide by the terms and conditions herein. You agree that you are solely responsible for your conduct while using this product, and for any consequences thereof. You agree to use this product only for proper purposes and all applicable laws, rules, and regulations, and all terms, precautions, practices, policies, and guidelines. Our company has made and may make available. Our company accepts no liability for damage, injury or, any legal responsibility incurred directly or indirectly from the use of the product. Our company will not provide any service for any product obtained from abnormal channels.

Safety Instruction

1. This device is not recommended to directly massage the skin to avoid harming the skin. Customer should massage the body after wearing dry clothes or placing a clean and dry towel on massage part, gently press and move the device.
2. High-speed and high-pressure massage can only be used to soft tissues without feeling pain or maladjustments. Do not massage head and any skeletal sites with high speed and high pressure.
3. Massage the same part of the body with high speed and high pressure for a long time will cause bruises. Please stop using immediately after feeling pain or maladjustments.

4. Please keep your fingers, hair or, any other body hair away from the extendable part of the device to prevent pinching and clipping.
5. Do not discard or misuse this device.
6. Do not modify this device by yourself without permission.
7. Do not operate or charge the device when nobody was around.
8. This device is for adults only.
9. People with contraindicated symptoms should use this device with doctor's permission.

Notice

1. Make sure the massage head is not blocked by external force when the device is power on. If it is blocked, it will pause for a while and try to restart.
2. DO NOT contact the product with water or other liquid for the Percussion Massage Device is not waterproof.
3. DO NOT disassemble the product except marked detachable. It need send to after-sales or authorized service center to fix it if you accidentally disassemble it and cause abnormal work. The relevant costs are borne by users.
4. Prolonged continuous operation may cause the motor surface temperature to rise, please operate carefully.
5. DO NOT drop or strike the product. If the product is abnormal, contact After-sales support.

Storage and Maintenance

1. Keep the product out of the reach of children and pets.
2. DO NOT leave the product near heat sources such as a furnace or heater. DO NOT leave the product inside of a vehicle on hot days.
3. Please storage the product in a dry environment.
4. DO NOT overcharge or overuse the battery, otherwise, it will cause damage to the battery core.
5. Never use the product when the temperature is too high or too low.

This document is subject to change without notice.

If you have any questions about this document, please contact us by sending a message to service@kica-tech.com.

Manufactured by: Guilin Feiyu Technology Incorporated Company

Website: www.kica-tech.com | Tel: +86-773-2320865 | E-mail: service@kica-tech.com