
Foot Massager

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
V 1.0

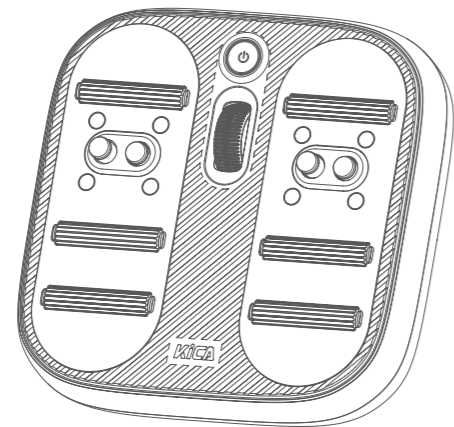EN Instructions


1. Overview

- [1] Power button/ Indicator
- [2] Speed control wheel
- [3] Massage head
- [4] Massage stick
- [5] Magnetite

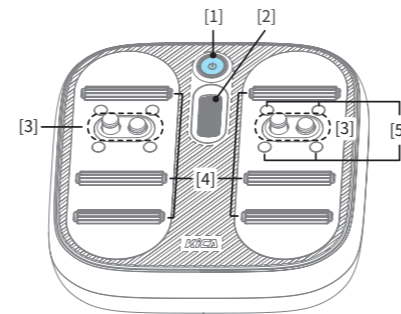2. Power on/Power off/Standby

Long press power button to power on/off.
 Click power button once after power on to enter standby mode, click again to wake up.

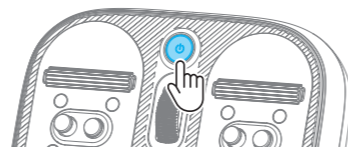3. Speed control

Roll the speed control wheel to adjust revolution speed of massage heads after power on. Roll upwards to improve speed, and roll downwards to reduce speed.

The massager default to be in high-frequency percussion mode after power on, and it enter gent clapping mode after the speed drop below 450 rpm.

Mode	Rated revolution speed (rpm)
High-frequency percussion mode	450~2200
Gent clapping mode	50~450

4. Locked motor protection

If the rotation is blocked for 3 seconds, it will enter standby mode automatically after a beep sound along with the red light flashing.

- ⚠ Do not press the massage head firmly to block the motor rotation while using, otherwise it will reduce the service life of the product! !
- When the battery is low, press the massage head firmly to block the motor rotation will easily cause the product to automatically power off.

5. Intelligent Protection

In order to prevent injury caused by massaging for a long time, it will enter standby mode after working for 20 minutes continuously.

If there is no operation for 3 minutes after enter standby mode, it will power off automatically.


6. Indicator

Indicator when not charging	Description
Constant green light	Working properly
Breathing green light	Standby
Green light flash once	Reach the maximum/Minimum revolution speed
Red light flashes slowly	Low power
Red light flashes quickly	Low power, will auto power off/Overheating protection/Motor locked

Indicator while charging	Description
Constant red light	Charging while power off
Constant yellow light	Charging while power on
Breathing yellow light	Charging while standby
Red light flash once	Reach the maximum/Minimum revolution speed
Constant green light	Charging completed
Red light flashes quickly	Low power, will auto power off/Overheating protection/Motor locked

7. Charging

Charging with USB-C cable.

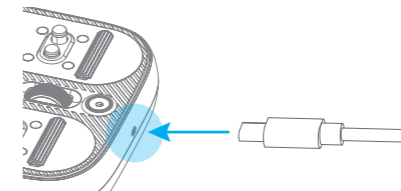8. Specifications

Product name: KiCA Foot Massager
 Product model: KiCAMG01
 Operating voltage range: 14V-16.8V
 Operating temperature range: 0°C ~+40°C
 Storage temperature range: -20°C ~+60°C
 Battery capacity: 2200mAh
 Rated input: 5V-15V ⚡ , 1A-2A
 Charging time: About 3 hours
 Weight: About 1890g
 Size: 300×300×83.7mm
 Max battery life without load: > 18 hours

List

Foot Massager×1, USB-C Cable ×1, Manual×1

Contraindications

1. Cancers
2. Metal bone anastomosis (do not use the device directly in the fracture area)
3. Osteoporosis
4. Vein thrombosis and pulsation
5. Hypertension
6. Pregnant woman
7. Varicose veins
8. Acute inflammation
9. Do not massage the bony area
10. Do not massage surgical site
11. Do not massage the same area for over 5 minutes

Safety Instruction

1. Keep the massage area dry and clean. It is recommended to use with shoe covers or socks, and you need to press and move gently.
2. In the case of not producing pain and discomfort, only the soft tissue parts of the foot can be used at high speed and high pressure, and it is not allowed to use high speed and high pressure on any skeletal parts.
3. Massage the same part of the body with high speed and high pressure for a long time will.
4. Please keep your fingers, hair or, any other body hair away from the extendable part of the device to prevent pinching and clipping.
5. Do not discard or misuse this device.
6. Do not modify this device by yourself without permission.
7. Do not operate or charge the device when nobody was around.
8. This device is for adults only.
9. People with contraindicated symptoms should use this device with doctor's permission.

Notice

1. Make sure the massage head is not blocked by external force when the device is power on. If it is blocked, it will pause for a while and try to restart.
2. DO NOT contact the product with water or other liquid for the massage device is not waterproof.
3. DO NOT disassemble the product except marked detachable. It need send to after-sales or authorized service center to fix it if you accidentally disassemble it and cause abnormal work. The relevant costs are borne by users.
4. Prolonged continuous operation may cause the motor surface temperature to rise, please operate carefully.
5. DO NOT drop or strike the product. If the product is abnormal, contact After-sales support.

Storage and Maintenance

1. Keep the product out of the reach of children and pets.
2. DO NOT leave the product near heat sources such as a furnace or heater. DO NOT leave the product inside of a vehicle on hot days.
3. Please storage the product in a dry environment.
4. DO NOT overcharge or overuse the battery, otherwise, it will cause damage to the battery core.
5. Never use the product when the temperature is too high or too low.

Precautions when discarding the product

The lithium ion battery used by the product is a recyclable and renewable valuable resource. When discarding it, please have it handled by a professional agency that can recycle chargeable batteries. Before discarding the product, power it off and take the lithium ion battery out from it following the steps below. The battery is to be disposed of safely.

1. Use a phillips screwdriver or other tools to open or pry off the back cover of the product.
2. Pull out the batter terminals from the PCBA.
3. Take out the battery.

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This document is subject to change without notice.
说明书内容如有更新，恕不另行通知。

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Warranty:

1. Within one year from the date of selling, the product is malfunctioning under normal due to non-artificial reasons.
2. The malfunction of the product is not caused by artificial reasons such as unauthorized disassemble or conversion or addition.
3. The buyer is able to provide the certificate of maintenance service: the warranty card, the legitimate receipts, invoices, or screenshot of purchasing.

The followings cases are not included under the warranty:

1. Unable to provide legitimate receipt and warranty card with buyer's info.
2. The damage caused by human or irresistible factor.

More details about after-sales policy, please refer to after-sales page on website:
<https://www.kica-tech.com/en/support>

* Our company reserves the right of final interpretation on the above-mentioned after sales terms and limitations.

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