

KiCA[®]


MINI-C

MESSAGE DEVICE

Instructions **V1.0**

Overview

1. Flat Head (Non-removable)
2. Battery Indicator
3. Charging Port
4. Power Button (Selecting gears/Power switch)
5. Gear Indicator (Breathing light)

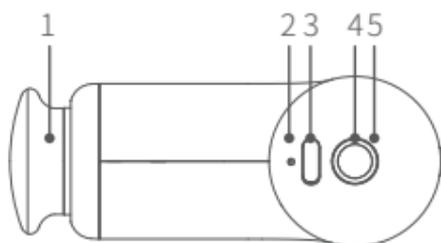

Flat Head

Massage flat major muscle groups on the back.



Spherical Head

Relax soft tissues and sensitive muscles.



Cone Head

Percuss deep tissues like meridians and acupoint, etc.




U-shaped Head

Massage neck (not cervical vertebrae) and tendon.

Start to Use

1. Massage Head Combination


The device can be used directly with the built-in massage head (flat head) on it, or combine with another one depending on your needs. See complete combination as below.

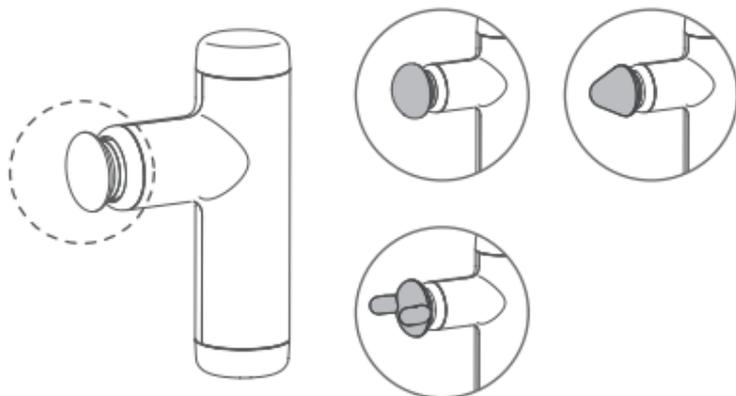

2. Power ON/OFF

Power on: Long press the power button.

Power off: Long press power button or single tap it to switch modes until the product power off.

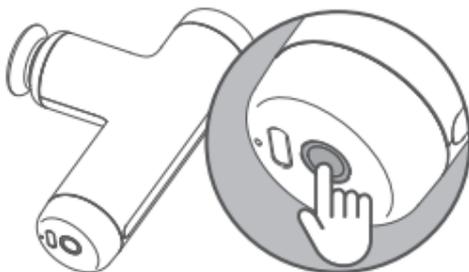

3. Selecting Gears

Single tap power button to select gears after powering on.

Gears	Rated rotate speed (rpm)	Function
First gear	1500	Muscle arousal
Second gear	2400	Relax
Third gear	3200	Deeply massage
Fourth gear	650	Gently massage
Power off	0	

Gears	Gear indicator (breathing light)
First gear	Blue breathing light
Second gear	Green breathing light
Third gear	Red breathing light
Fourth gear	Yellow breathing light

4. Intelligent Protection

This product has intelligent protection function. In order to prevent muscle injury caused by massaging for a long time, a beeping sound will be issued and gear indicator turn off (battery indicator remain) after working for 10 minutes continuously, and it can be restarted by pressing the function button. If there is no operation within 5 seconds, it will power off automatically.

5. Locked Motor Protection

This product has Locked Motor Protection function. If the rotating is blocked for 3 seconds, the product will power off automatically.

 Do not press the massage head firmly to block the motor rotation while using, otherwise it will reduce the service life of the product!

6. Battery Indicator

Indicator while using	Description
Constant green light	Fully charged
Green light flashes slowly	Medium level of battery
Green light flashes quickly	Low power

Indicator while charging	Description
Constant red light	Charging
Constant green light	Charge completed

7. Charging

⚠ Do not use it while charging.



8. Specifications

Product Name: KiCA Massage Device

Product Model: MINI-C

Operating Voltage Range: 7.1V~8.4V

Operating Temperature Range: 0°C~+40°C

Charging Temperature Range: 0°C~+40°C

Battery Capacity: 1300mAh

Charging Time: 2h

Weight: 213g

Size: 88.2×35×128mm

Theoretical Battery Life: 10 Hours (Test in lab environment under the fourth gear with flat head)

List

- | | |
|-----------------------------------|----|
| ① Main Body (including flat head) | ×1 |
| ② Spherical Head | ×1 |
| ③ Cone Head | ×1 |
| ④ U-shaped Head | ×1 |
| ⑤ USB-C Cable | ×1 |
| ⑥ Manual | ×1 |



Contraindications

1. Cancers.
2. Metal bone anastomosis (do not use the device directly in the fracture area).
3. Osteoporosis.
4. Vein thrombosis and pulsation.
5. Hypertension.
6. Pregnant woman.
7. Varicose veins.
8. Acute inflammation.
9. Do not massage the bony area.
10. Do not massage surgical site.
11. Do not massage the cervical vertebrae and spine.
12. Kidney stones (do not use the device near the kidneys).
13. Do not massage one episode for more than 5 minutes.

Disclaimer

Thanks for using the product. The information in this document affects your safety, legal rights, and responsibilities. Read the entire document carefully to ensure proper configuration before using. Failure to read and follow instructions and warnings in this document may result in serious injury to yourself or others, or damage your products or other objects in the vicinity.

By using this product, you hereby signify that you have read this

disclaimer and warning carefully and that you understand and agree to abide by the terms and conditions herein. You agree that you are solely responsible for your conduct while using this product, and for any consequences thereof. You agree to use this product only for proper purposes and all applicable laws, rules, and regulations, and all terms, precautions, practices, policies, and guidelines. Our company has made and may make available. Our company accepts no liability for damage, injury or, any legal responsibility incurred directly or indirectly from the use of the product.

Our company will not provide any service for any product obtained from abnormal channels.

Disclaimer

1. It is recommend to massage your body after wearing dry clothes or placing a clean and dry towel on massage part, then gently press and move the device.
2. High-speed and high-pressure massage can only be used to soft tissues without feeling pain or maladjustments. Do not massage head and any skeletal sites with high speed and high pressure.
3. Massage the same part of the body with high speed and high pressure for a long time will.
4. Please keep your fingers, hair or, any other body hair away from the extendable part of the device to prevent pinching and clipping.
5. Do not discard or misuse this device.
6. Do not modify this device by yourself without permission.
7. Do not operate or charge the device when nobody was around.
8. This device is for adults only.
9. People with contraindicated symptoms should use this device with doctor's permission.

Notice

1. Make sure the massage head is not blocked by external force when the device is power on. If it is blocked, it will pause for a while and try to restart.
2. DO NOT contact the product with water or other liquid for the massage device is not waterproof.
3. DO NOT disassemble the product except marked detachable. It need send to after-sales or authorized service center to fix it if you accidentally disassemble it and cause abnormal work. The relevant costs are borne by users.
4. Prolonged continuous operation may cause the motor surface temperature to rise, please operate carefully.
5. DO NOT drop or strike the product. If the product is abnormal, contact After-sales support.

Storage and Maintenance

1. Keep the product out of the reach of children and pets.
2. DO NOT leave the product near heat sources such as a furnace or heater. DO NOT leave the product inside of a vehicle on hot days.
3. Please storage the product in a dry environment.
4. DO NOT overcharge or overuse the battery, otherwise, it will cause damage to the battery core.
5. Never use the product when the temperature is too high or too low.

Precautions when discarding the product

The lithium ion battery used by the product is a recyclable and renewable valuable resource.

When discarding it, please have it handled by a professional agency that can recycle chargeable batteries. Before discarding the product, power it off and take the lithium ion battery out from it following the steps below. The battery is to be disposed of safely.

1. Disassemble the top housing with help of a professional tool.
2. Pull out the batter terminals from the PCBA.
3. Use a professional cutter to cut off the bottom housing and then take out the battery.

Manufactured by: Guilin Feiyu Technology Incorporated Company
Address: Feiyu Intelligent Industrial Park, D-14-1b, Information Industry Park, National High-tech Zone, ChaoYang Road, Qi Xing District, Guilin, Guangxi, China
Website: www.kica-tech.com
Tel: 0773-2320865
E-mail: service@feiyu-tech.com

The logo for KICA, featuring the word "KICA" in a bold, teal, sans-serif font. A registered trademark symbol (®) is positioned to the upper right of the letter "A". The logo is centered within a large, light teal outline of a stylized letter "A" that serves as a background element.

KICA[®]

www.kica-tech.com