


# ***KiCAMINI***



## **Massage Device**

Instructions (V2.1)

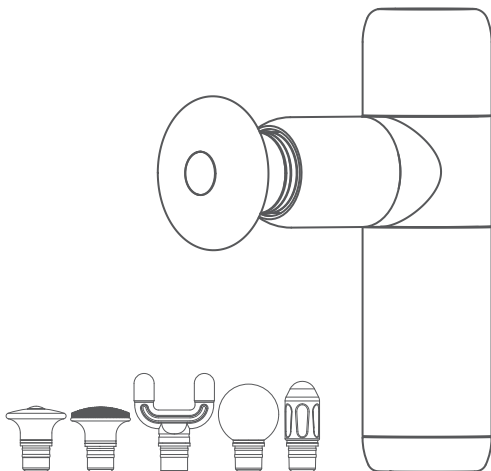

## Overview

1. Massage Head (Replaceable)
2. Battery Indicator
3. Charging Port
4. Power Button  
(Selecting gears/Power switch)
5. Gear Indicator (Breathing light)

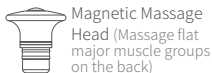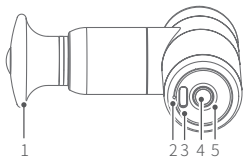


Magnetic Massage Head (Massage flat major muscle groups on the back)



Gentle Head (Promote the absorption of essential oil)



Spherical Head (Relax soft tissues and sensitive muscles)

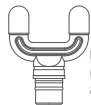

U-shaped Head (Massage neck and tendon)




Cone Head (Percuss deep tissues like meridians and acupoint, etc.)

## Start to Use

### 1. Install/ Replace Massage Head

Install a massage head into the hole in front of the device according to your need.

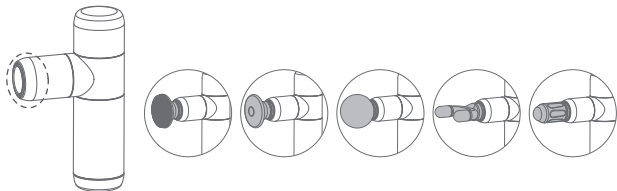


### 2. Power ON

Long press the Power button.

### 3. Power OFF

Single tap Power button to switch to standby mode or long press Power button.

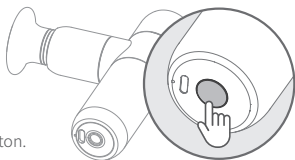

## 4. Selecting Gears

Single tap power button to select gears after powering on.

Gears	Gear Indicator (breathing light)	Rated Rotate Speed (rpm)	Function
First gear	Blue breathing light	1500	Muscle arousal
Second gear	Green breathing light	2400	Relax
Third gear	Red breathing light	3200	Deeply massage
Fourth gear	Yellow breathing light	650	Gently massage
Standby mode (Auto power off after 2 seconds)		0	

## 5. Intelligent Protection

This product has intelligent protection function. In order to prevent muscle injury caused by massaging for a long time, a beeping sound will be issued and gear indicator turn off ( battery indicator remain ) after working for 10 minutes continuously, and it can be restarted by pressing the function button. If there is no operation within 5 seconds, it will power off automatically.

## 6. Locked Motor Protection

This product has Locked Motor Protection function. If the rotating is blocked for 3 seconds, the product will enter standby mode automatically. Users can restart it by pressing the function button. If there is no operation in 10 seconds, the product will power off automatically.

## 7. Gentle head introduction

Gentle head is used to gently massage sensitive part of body with gentle mode. Gentle head is suitable for oil massage with the silicone brush design which can promote the absorption of essential oil.




## 8. Battery Indicator

Indicator while using	Description
Constant green light	Fully charged
Green light flashes slowly	Medium level of battery
Green light flashes quickly	Low power

Indicator while charging	Description
Constant red light	Charging
Constant green light	Charge completed

## 9. Charging

!Do not use it while charging

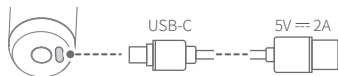

## 10. Specifications

Product name: KiCA Massage Device

Model: MINI

Operating Voltage Range: 7.1V~8.4V


Battery Capacity: 1300mAh

Charging Time: >2h

Weight (Reference value): 248g

Size: 77×35×125.6mm

Theoretical Battery Life: About 8 Hours (Test in lab environment under the fourth gear)

 Do not press the massage head firmly to block the motor rotation while using, otherwise it will reduce the service life of the product!

### List

- |                         |    |
|-------------------------|----|
| ① Main Body             | x1 |
| ② Magnetic Massage Head | x1 |
| ③ Gentle head           | x1 |
| ④ Cone Head             | x1 |
| ⑤ Spherical Head        | x1 |
| ⑥ U-shaped Head         | x1 |
| ⑦ USB-C Cable           | x1 |
| ⑧ Storage Bag           | x1 |
| ⑨ Manual                | x1 |



## Disclaimer

Thanks for using the product. The information in this document affects your safety, legal rights, and responsibilities. Read the entire document carefully to ensure proper configuration before using. Failure to read and follow instructions and warnings in this document may result in serious injury to yourself or others, or damage your products or other objects in the vicinity.

By using this product, you hereby signify that you have read this disclaimer and warning carefully and that you understand and agree to abide by the terms and conditions herein. You agree that you are solely responsible for your conduct while using this product, and for any consequences thereof. You agree to use this product only for proper purposes and all applicable laws, rules, and regulations, and all terms, precautions, practices, policies, and guidelines. Our company has made and may make available. Our company accepts no liability for damage, injury or, any legal responsibility incurred directly or indirectly from the use of the product. Our company will not provide any service for any product obtained from abnormal channels.

## Safety Instruction

1. This device is not recommended to directly massage the skin except gentle head to avoid to do harm to the skin. Customer should massage body after wearing dry clothes or placing a clean and dry towel on massage part, gently press and move the device.
2. High-speed and high-pressure massage can only be used to soft tissues without feeling pain or maladjustments. Do not massage head and any skeletal sites with high speed and high pressure.
3. Massage the same part of the body with high speed and high pressure for a long time will.
4. Please keep your fingers, hair or, any other body hair away from the extendable part of the device to prevent pinching and clipping.
5. Do not discard or misuse this device.
6. Do not modify this device by yourself without permission.
7. Do not operate or charge the device when nobody was around.
8. This device is for adults only.
9. People with contraindicated symptoms should use this device with doctor's permission.

## Notice

1. Make sure the massage head is not blocked by external force when the device is power on. If it is blocked, it will pause for a while and try to restart.
2. DO NOT contact the product with water or other liquid for the Massage Device is not waterproof.
3. DO NOT disassemble the product except marked detachable. It need send to after-sales or authorized service center to fix it if you accidentally disassemble it and cause abnormal work. The relevant costs are borne by users.
4. Prolonged continuous operation may cause the motor surface temperature to rise, please operate carefully.
5. DO NOT drop or strike the product. If the product is abnormal, contact After-sales support.

## Storage and Maintenance

1. Keep the product out of the reach of children and pets.
2. DO NOT leave the product near heat sources such as a furnace or heater. DO NOT leave the product inside of a vehicle on hot days.
3. Please storage the product in a dry environment.
4. DO NOT overcharge or overuse the battery, otherwise, it will cause damage to the battery core.
5. Never use the product when the temperature is too high or too low.

## Precautions when discarding the product

The lithium ion battery used by the product is a recyclable and renewable valuable resource. When discarding it, please have it handled by a professional agency that can recycle chargeable batteries. Before discarding the product, power it off and take the lithium ion battery out from it following the steps below. The battery is to be disposed of safely.

1. Disassemble the top housing with help of a professional tool.
2. Pull out the batter terminals from the PCBA.
3. Use a professional cutter to cut off the bottom housing and then take out the battery.

## Contraindications

- (1) Cancers.
- (2) Metal bone anastomosis (do not use the device directly in the fracture area).
- (3) Osteoporosis.
- (4) Vein thrombosis and pulsation.
- (5) Hypertension.
- (6) Pregnant woman.
- (7) Kidney stones (do not use the device near the kidneys).
- (8) Varicose veins.
- (9) Acute inflammation.
- (10) Do not massage one episode for more than 5 minutes.
- (11) Do not massage the bony area.
- (12) Do not massage surgical site .
- (13) Do not massage the cervical vertebrae and spine.

## Find More Features



Website



Youtube



Facebook

This document is subject to change without notice. If you have any questions about this document, please contact us by the following ways.

Manufactured by: Guilin Feiyu Technology Incorporated Company

Website: [www.kica-tech.com](http://www.kica-tech.com)

Tel: +86-773-2320865

E-mail: [service@kica-tech.com](mailto:service@kica-tech.com)