


KiCA Pro

Massage Device Instructions

V1.0

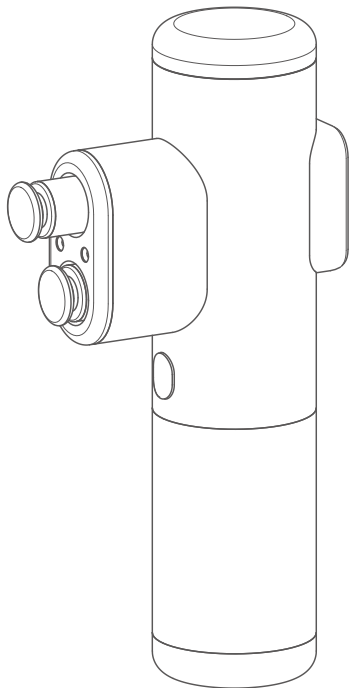


Contents

Overview	1
Start to Use	1
1. Massage Head Combination	1
2. Power ON / Start running / Power Off	2
3. Touch screen	3
3.1 Slide from right to left in home page	5
3.2 Slide right in home page	6
3.3 Slide down in home page	7
3.4 Slide up in home page	7
4. Locked Motor Protection	8
5. Intelligent Protection	8
6. Charging	8
7. Specifications	9

Overview

KiCA Pro is a massage device which developed by Guilin Feiyu Technology Incorporated Company. It comes with double massage heads and touch screen. KiCA Pro provide comfortable and convenient massage experience for users.

- [1] Flat head (Non-removable)
- [2] Trigger button
- [3] Touch screen
- [4] Power button
- [5] Charging port
- [6] Battery indicator

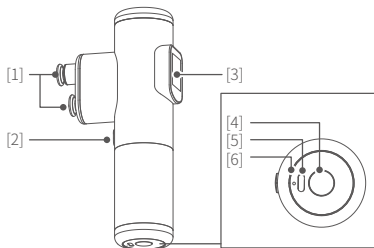

Flat Head:

Massage flat major muscle groups on the back.



Gentle Head:

Promote the absorption of essential oil.



Spherical Head:

Relax soft tissues and sensitive muscles.




U-shaped Head:

Massage neck (not cervical vertebrae) and tendon.

Start to Use

1. Massage Head Combination

The device can be used directly with the built-in massage head (two flat head), on it, or combine with another one depending on your needs. See complete combination as below.

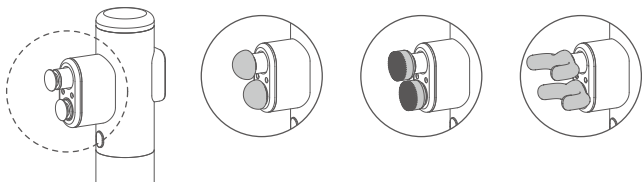


2. Power ON / Start running / Power Off

Power ON: Long press the power button until hear the tone.




Start running:

Method one: Click trigger button to start running, click again to stop.

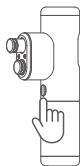


* Default to start running via "click" trigger button, it can also be customized as "press" in setting page. (3.2)

Method two: Operate on the touch screen to select the modes and gears.

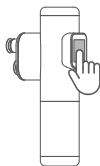



Do not touch the middle and the back of the two massage heads while the product is running to avoid injuries.

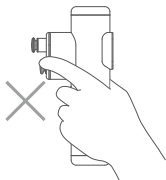

Power Off: Long press the power button until hear the tone to power off. If there is no operation within 10 minutes, it will power off automatically.

3. Touch screen

Operate on the touch screen to select the modes and gears.

- [1] Battery level
- [2] Current rotation speed
- [3] Massage time for this time
- [4] Total massage time for current mode
- [5] Current gear (1-10)
- [6] Current mode
- [7] Current pressure



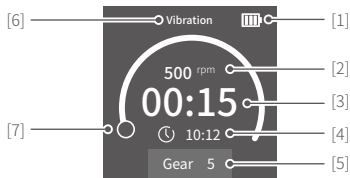


Diagram:



Basic message modes

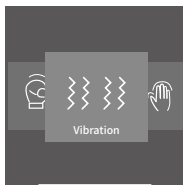

Setting

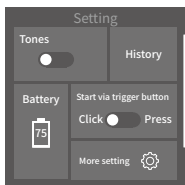


Package choosing

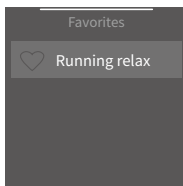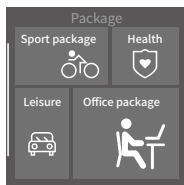


Favorites

Tap "Gear" icon to enter gear setting, slide to switch among 1-10 gears.

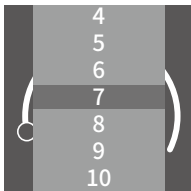

3.1 Slide from right to left in home page


This product preinstall several massage package, single tap corresponding icon to enter.

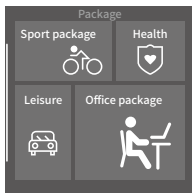

Package	Option	Introduction
Sport package	9	Running relax / Basketball / Soccer relax / Swimming relax / Relax leg / Relax chest / Relax waist / Relax hunkers / Relax back
Health	4	Cervical vertebrae / Lumbar vertebrae / Around shoulder / Knee-joint
Leisure	4	Driver / Office worker / Fat reduction / Consolidate kidney
Office package	3	Sedentary work / Longtime standing / Refreshing / Sleep

Take "Sport package-Running relax" as an example:

Tap ♡ icon to add corresponding mode in Favorites (Slide up in home page) (3.4)。

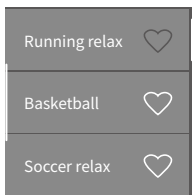

Tap "Running relax" to enter

[1] The sequence number for current message part/Total number of message parts

[2] Current message part name


[3] Massage time

[4] Diagram for current message part

[5] Last message part

[6] Current message state


[7] Next message part

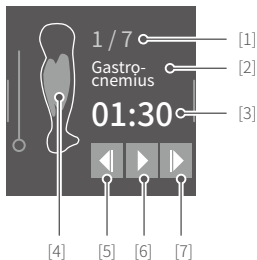

Operation:

Slide right——Exit the current page.

Slide left——Single tap to choose massage mode, slide to switch gears. Slide right in the middle of screen to exit after complete the setting.

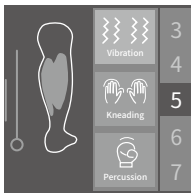

Usage:

· Tap [5] or [7] on touch screen to leaf through message part.

· Tap [6] to start (or click trigger button).

· Massage your body as [4] shown.

· The device will auto stop after finish the massage time which shown in [3]. Or click trigger button to stop midway.




When massaging neck and shoulders, make sure your hair will not be entangled into the piston.

3.2 Slide right in home page


Set tones, start via trigger button and more settings in this page.

You can also check history and battery level in this page.

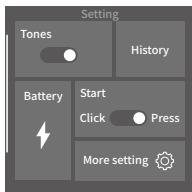

3.3 Slide down in home page


Basic massage modes:

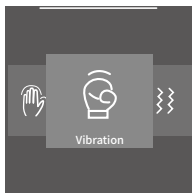

Mode	Rated rotate speed (rpm)	Gears
Vibration mode	500-2600	1-10
Kneading mode	100-900	1-10
Percussion mode	50-450	1-10

3.4 Slide up in home page

Show the modes which you have added to Favorites.

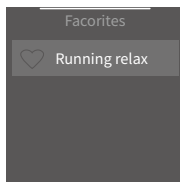

Side left and tap "Del" to remove from the Favorites.



4. Locked Motor Protection



Do not press the massage head firmly to block the motor rotation while using, otherwise it will reduce the service life of the product!

This product has locked motor protection function. If the rotation is blocked for 3 seconds, it will stop running automatically after a beep sound. Single tap the trigger button to start again.

5. Intelligent Protection

This product has intelligent protection function. After working for 10 minutes continuously, a beep sound will be issued, and the product will stop running automatically.

Single tap the trigger button to start again.

6. Charging




Do not use it while charging.


The battery indicator is red while charging, and it will turn green after charging completed. Screen shows the charging icon.

When the product is on, it will return to the interface before charging, after unplug the power supply while charging.




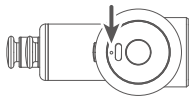

Battery indicator

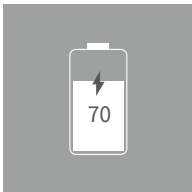

Charging

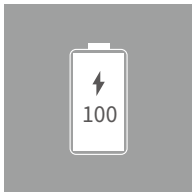

Charging completed

7. Specifications

Product Name: KiCA Pro Massage Device

Product Model: A3D

Operating Voltage Range: 10.8V~12.6V

Operating Temperature Range: 0°C ~+40°C

Charging Temperature Range: 0°C ~+40°C

Battery Capacity: 2200mAh

Charging Time: about 4.5Hours

Weight: about 630g

Size: 97.2×46×188.3mm

Theoretical Battery Life: About 20Hours (Test in lab environment under the 1st gear with flat head)

List

- | | |
|-----------------------------------|----|
| ① Main Body (including flat head) | ×1 |
| ② Storage Bag | ×1 |
| ③ Spherical Head | ×2 |
| ④ Gentle Head | ×2 |
| ⑤ U-shaped Head | ×2 |
| ⑥ USB-C Cable | ×1 |
| ⑦ Manual | ×1 |



Contraindications

1. Cancers
2. Metal bone anastomosis (do not use the device directly in the fracture area)
3. Osteoporosis
4. Vein thrombosis and pulsation
5. Hypertension
6. Pregnant woman
7. Varicose veins
8. Acute inflammation
9. Do not massage the bony area
10. Do not massage surgical site
11. Do not massage the cervical vertebrae and spine
12. Kidney stones (do not use the device near the kidneys)
13. Do not massage one episode for more than 5 minutes

Safety Instruction

1. It is recommend to massage your body after wearing dry clothes or placing a clean and dry towel on massage part, then gently press and move the device.
2. High-speed and high-pressure massage can only be used to soft tissues without feeling pain or maladjustments. Do not massage head and any skeletal sites with high speed and high pressure.
3. Massage the same part of the body with high speed and high pressure for a long time will.
4. Please keep your fingers, hair or, any other body hair away from the extendable part of the device to prevent pinching and clipping.
5. Do not discard or misuse this device.
6. Do not modify this device by yourself without permission.
7. Do not operate or charge the device when nobody was around.
8. This device is for adults only.
9. People with contraindicated symptoms should use this device with doctor's permission.

Notice

1. Make sure the massage head is not blocked by external force when the device is power on. If it is blocked, it will pause for a while and try to restart.
2. DO NOT contact the product with water or other liquid for the massage device is not waterproof.
3. DO NOT disassemble the product except marked detachable. It need send to after-sales or authorized service center to fix it if you accidentally disassemble it and cause abnormal work. The relevant costs are borne by users.
4. Prolonged continuous operation may cause the motor surface temperature to rise, please operate carefully.
5. DO NOT drop or strike the product. If the product is abnormal, contact After-sales support.

Storage and Maintenance

1. Keep the product out of the reach of children and pets.
2. DO NOT leave the product near heat sources such as a furnace or heater. DO NOT leave the product inside of a vehicle on hot days.
3. Please storage the product in a dry environment.
4. DO NOT overcharge or overuse the battery, otherwise, it will cause damage to the battery core.
5. Never use the product when the temperature is too high or too low.

Precautions when discarding the product

The lithium ion battery used by the product is a recyclable and renewable valuable resource.

When discarding it, please have it handled by a professional agency that can recycle chargeable batteries. Before discarding the product, power it off and take the lithium ion battery out from it following the steps below. The battery is to be disposed of safely.

1. Disassemble the top housing with help of a professional tool.
2. Pull out the batter terminals from the PCBA.
3. Use a professional cutter to cut off the bottom housing and then take out the battery.

Find More Features



Website



Youtube



Facebook



This document is subject to change without notice.

If you have any questions about this document, please contact us by the following ways.

Manufactured by: Guilin Feiyu Technology Incorporated Company

Website: www.kica-tech.com

Tel: +86 773-2320865

E-mail: service@kica-tech.com
